

# Anti-HA Mouse Monoclonal Antibody

Technical literature is available at: [www.mesgenbio.com](http://www.mesgenbio.com). E-mail MesGen Technical Services if you have questions on use of this system: [tech@mesgenbio.com](mailto:tech@mesgenbio.com)

**MESGEN**  
INNOVATION BIOTECHNOLOGY

Catalog Number : MAN2003

Lot Number : Refer to vial

Packaging Size : 100µL  200µL  500µL  1000µL

## Background

Recognizes HA-tagged proteins overexpressed in cells, including both amino & carboxy termini of targeted proteins in transfected mammalian cells.

## Overview

**Description** : Mouse monoclonal antibody to HA-Tag [6A1]

**Tested applications** : WB, ICC/IF, IP

**Formulation** : Liquid


## Properties

**Immunogen** : YPYDVPDYA conjugated to KLH.

**Clonality** : 6A1

**Isotype** : Mouse IgG2

## Application notes

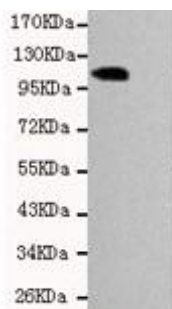

**Fig.** Western blot using of 293 cells transfected with HA-tagged vector(1) at 1:5000 and untransfected control (2).

Optimal working dilutions should be determined experimentally by the investigator. Suggested starting dilutions are as follows:  
WB (1:1000-1:10000)                      ICC/IF (1:100-1:500)                      IP (1:50-200)

## Storage instruction

Store at +4°C short term (1-2 weeks).

Aliquot and store at -20°C or -80°C

**Note:** The product listed herein is for research use only and is not intended for use in human or clinical diagnosis. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license. We cannot be responsible for patent infringements or other violations that may occur with the use of this product.

## Distributor in China ( Mainland )

**Company** : Shanghai hongsheng Biotechnology Co.,Ltd

**Address** : Room 2-205, Xingzhan building, No.185, Liuying Road, Zhabei District, Shanghai

**Tel** : 86-21-56620378    **Fax** : 86-21-56620378

**Email** : [sales@hongshengbio.com](mailto:sales@hongshengbio.com)

**Website** : [www.hongshengbio.com](http://www.hongshengbio.com)